

# AccuCore<sup>®</sup>

## Single Action Biopsy Needles



Fully Exposed Notch  
For easy sample removal



Matching Coaxial Needle  
With HiLiter<sup>®</sup> ultrasound enhancement



Trocar Tip (14 gauge Only)  
For easy insertion



...helping to reduce the uncertainty in diagnosis.




## Stronger Spring

For increased cutting effectiveness the AccuCore® is constructed with a stronger spring to improve specimen acquisition.



## Trocar Point

The 14 gauge AccuCore® features a Trocar point for ease of insertion into dense tissue. Other gauge sizes feature a standard bevel tip.

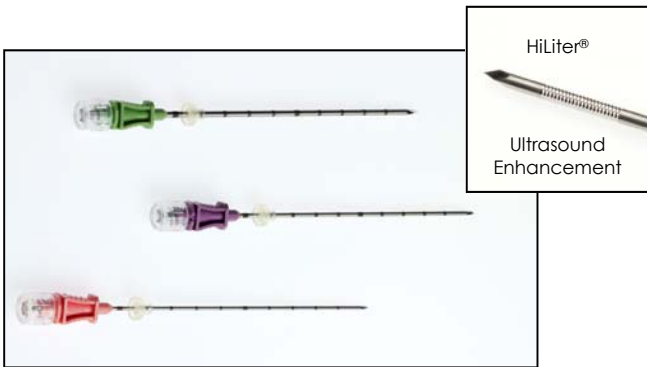

## Matching Coaxial Needles

Matching coaxial needles for each AccuCore® size are available and are color-coded by size. Coaxial needles feature the patented HiLiter® enhancement which significantly increases visibility under ultrasound.

## Additional Features

- Lightweight and compact—Easy to Maneuver
- Sharp cutting cannula for intact core sample
- Specimen notch is fully exposed when stylet is advanced for ease of sample removal
- Ultrasound enhanced tip for easier visualization under ultrasound guidance
- Depth marks on cannula for precise positioning



Ordering Information					
AccuCore® Biopsy Device					
Code	Gauge	Length (CM)	Length (In)	Tip	Coaxial w/HiLiter®
640914	14	9	3 ½"	Trocar	601409
641514	14	15	5 ⅞"	Trocar	601415
620916	16	9	3 ½"	Standard	601609
621516	16	15	5 ⅞"	Standard	601615
620918	18	9	3 ½"	Standard	601809
621518	18	15	5 ⅞"	Standard	601815
620920	20	9	3 ½"	Standard	602009
621520	20	15	5 ⅞"	Standard	602015

For More Information Call Toll Free:  
**1-800-558-4647**



4375 Donker Ct. SE  
Grand Rapids, MI 49512

Toll Free: 800-558-4647  
Tel: 616-301-7800  
Fax: 616-301-7799  
Email: [inrad-inc@inrad-inc.com](mailto:inrad-inc@inrad-inc.com)  
[www.inrad-inc.com](http://www.inrad-inc.com)

Complete information on INRAD devices, including indications, contraindications and cautions, is contained in the product package labeling and should be consulted before use.

Restricted Device: Federal law (USA) restricts the sale, distribution or use of this device to, by or on the lawful order of a physician. This product is designed as a disposable product for single patient use only.

©INRAD is a registered trademark of INRAD, Inc.  
©AccuCore is a registered trademark of INRAD, Inc.  
©HiLiter is a registered trademark of MEDTRONIC, Inc.  
© Copyright 2010 INRAD, Inc.  
06369D