

Precise Core™

Biopsy Device

*Designed for
precision biopsy
applications*



**No Forward
Throw**



Pre-Fire

Armed


Fired



**Arm and Fire
With One Hand**

INRAD®

...helping to reduce the uncertainty in diagnosis.

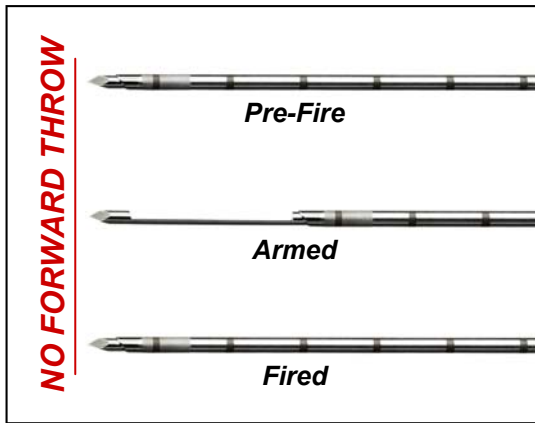


No Forward Throw

The PreciseCore™ needle does not advance when fired. Once positioned in the target area the device is armed, which draws back the outer cannula while the stylet remains stationary. This design offers precision and reduces the threat of damage to deeper tissue.



Trocar Point

The PreciseCore™ features a Trocar point for ease of insertion into dense tissue

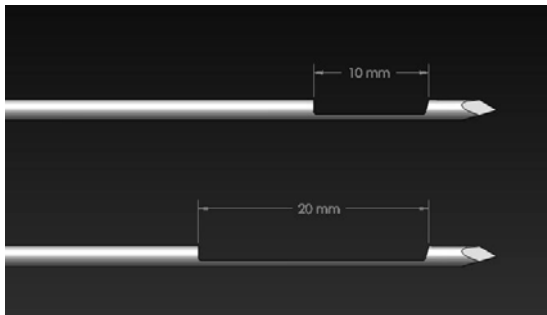

Two Sample Notch Options

For flexibility, the PreciseCore™ offers a choice of either a 10mm or 20mm notch which allows for greater precision and control depending on the procedure.

Additional Features

- Ergonomic design allows for ease of use along with comfort and control
- Easy to grip exterior and top-positioned arming button
- Forward-mounted trigger
- Ultrasound enhanced tip for increased visualization under ultrasound guidance
- Depth marks etched on cannula for precise positioning
- Color coded for easy gauge size identification



| Ordering Information | | | | |
|----------------------------|-------|-------------|-------------|--------------|
| PreciseCore™ Biopsy Device | | | | |
| Code | Gauge | Length (CM) | Length (In) | Notch Length |
| 670514 | 14 | 5 | 2" | 10mm |
| 671014 | 14 | 10 | 4" | 10mm |
| 670518 | 18 | 5 | 2" | 10mm |
| 671018 | 18 | 10 | 4" | 10mm |
| 670520 | 20 | 5 | 2" | 10mm |
| 671020 | 20 | 10 | 4" | 10mm |
| 680514 | 14 | 5 | 2" | 20mm |
| 681014 | 14 | 10 | 4" | 20mm |
| 680518 | 18 | 5 | 2" | 20mm |
| 681018 | 18 | 10 | 4" | 20mm |
| 680520 | 20 | 5 | 2" | 20mm |
| 681020 | 20 | 10 | 4" | 20mm |

For More Information Call Toll Free:
1-800-558-4647



4375 Donker Ct. SE
Grand Rapids, MI 49512

Toll Free: 800-558-4647
Tel: 616-301-7800
Fax: 616-301-7799
Email: inrad-inc@inrad-inc.com
www.inrad-inc.com

Complete information on INRAD devices, including indications, contraindications and cautions, is contained in the product package labeling and should be consulted before use.

Restricted Device: Federal law (USA) restricts the sale, distribution or use of this device to, by or on the lawful order of a physician. This product is designed as a disposable product for single patient use only.

©INRAD is a registered trademark of INRAD, Inc.
™PreciseCore is a trademark of INRAD, Inc.
© Copyright 2009 INRAD, Inc.
01901B