

# Precise Core™

## Biopsy Device

*Designed for  
precision biopsy  
applications*



**No Forward  
Throw**



*Pre-Fire*

*Armed*


*Fired*



**Arm and Fire  
With One Hand**

**INRAD®**

*...helping to reduce the uncertainty in diagnosis.*

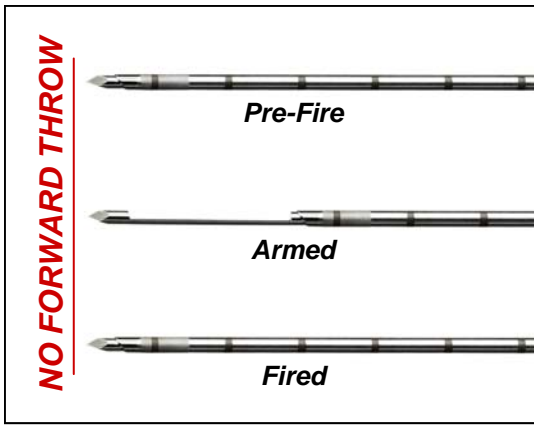


## No Forward Throw

The PreciseCore™ needle does not advance when fired. Once positioned in the target area the device is armed, which draws back the outer cannula while the stylet remains stationary. This design offers precision and reduces the threat of damage to deeper tissue.



## Trocar Point

The PreciseCore™ features a Trocar point for ease of insertion into dense tissue

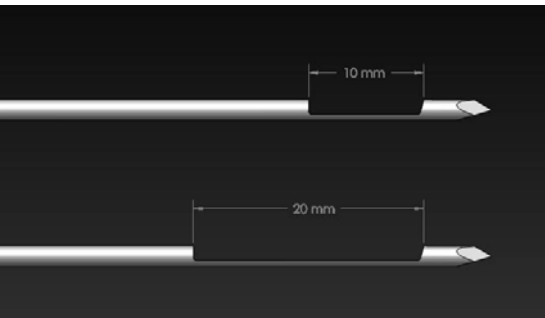

## Two Sample Notch Options

For flexibility, the PreciseCore™ offers a choice of either a 10mm or 20mm notch which allows for greater precision and control depending on the procedure.

## Additional Features

- Ergonomic design allows for ease of use along with comfort and control
- Easy to grip exterior and top-positioned arming button
- Forward-mounted trigger
- Ultrasound enhanced tip for increased visualization under ultrasound guidance
- Depth marks etched on cannula for precise positioning
- Color coded for easy gauge size identification



Ordering Information					
PreciseCore™ Biopsy Device					
Code	Gauge	Length (CM)	Length (In)	Notch Length	Coaxial w/ HiLiter® Tip
670514	14	5	2"	10mm	-
671014	14	10	4"	10mm	691014
670518	18	5	2"	10mm	-
671018	18	10	4"	10mm	691018
670520	20	5	2"	10mm	-
671020	20	10	4"	10mm	691020
680514	14	5	2"	20mm	-
681014	14	10	4"	20mm	601409
680518	18	5	2"	20mm	-
681018	18	10	4"	20mm	601809
681518	18	15	6"	20mm	691518
680520	20	5	2"	20mm	-
681020	20	10	4"	20mm	602009
681520	20	15	6"	20mm	691520

For More Information Call Toll Free:

**1-800-558-4647**



4375 Donker Ct. SE  
Grand Rapids, MI 49512

Toll Free: 800-558-4647

Tel: 616-301-7800

Fax: 616-301-7799

Email: [inrad-inc@inrad-inc.com](mailto:inrad-inc@inrad-inc.com)

[www.inrad-inc.com](http://www.inrad-inc.com)

Complete information on INRAD devices, including indications, contraindications and cautions, is contained in the product package labeling and should be consulted before use.

Restricted Device: Federal law (USA) restricts the sale, distribution or use of this device to, by or on the lawful order of a physician. This product is designed as a disposable product for single patient use only.

©INRAD is a registered trademark of INRAD, Inc.

™PreciseCore is a trademark of INRAD, Inc.

© Copyright 2009-2011 INRAD, Inc.

01901C